

Experiment 5

Simple Harmonic Motion

Reading:

Halliday, Resnick, Walker Chapter 14 as needed (e.g. definitions of oscillation variables)

No Homework, but Extra Credit is available:

1. Show by using Hooke's Law that a stationary mass m hanging from a spring with constant k (Fig. 1) stretches the spring to a new equilibrium position $x_0 = \frac{mg}{k}$.
2. Confirm by substitution that Eq. (7) indeed solves Eq. (6).
3. Show that when $\beta = 0$ (no damping), Eq. (7) reduces to Eq. (3).

1. Goals

1. Understand the properties of an oscillating system governed by Hooke's Law.
2. Study the effects of friction on an oscillating system, which leads to damping.
3. Use a non-linear least-squares fitting procedure to characterize an oscillator.
4. Learn the basics of using an oscilloscope.
5. Understand how to detect a power law relationship between variables.

2. Theoretical Introduction

2.1 Simple Harmonic Oscillation (SHO) Consider a system illustrated in the figure below. It consists of a mass m suspended from a spring with spring constant k .

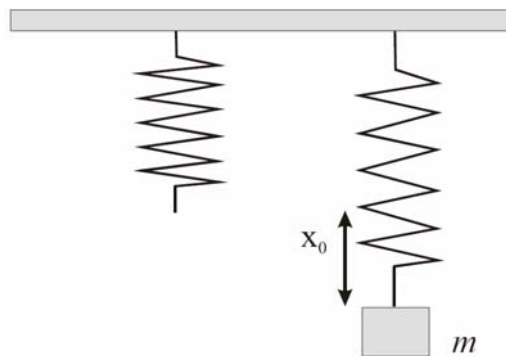


Fig. 1. A mass on a spring in the gravitational field of Earth

Hooke's law states that the force resisting the extension of the spring is proportional to the deviation of the spring from its equilibrium position. That is, $F = -kx$, where $x = 0$ is defined by the equilibrium position of the spring. If we gently add a mass to the spring, the spring will stretch to a new (lower) equilibrium position, $-x_0 = -mg/k$, where g is the gravitational constant. At this position, the vertical restoring force of the spring balances the weight. In what follows, we will take the origin of x at this new equilibrium position. In these final coordinates, at $x=0$ the gravitational force is canceled by the force due the spring, so that we can ignore the constant gravitational force in this coordinate system, and consider only additional Hooke's law forces beyond those required to cancel the weight.

If the mass is now displaced from its equilibrium position, the same equation $F = -kx$ still applies, where x is now that displacement. In other words, when such a displacement is made, a restoring

force again acts to return the mass to its equilibrium position. Upon release, the mass moves toward the equilibrium position, but its inertia causes it to “overshoot” this point. The motion then continues through the equilibrium position and beyond until the restoring force eventually stops the mass and pulls it back toward the equilibrium position. The motion then repeats itself back and forth through the position of equilibrium. Newton’s second law states that any unbalanced force results in an acceleration of the mass, proportional to the force. If we apply Newton’s second law to the motion of a mass m that is subject to Hooke’s law, we get

$$m \frac{d^2 x}{dt^2} = -kx, \quad (1)$$

which can be written as

$$\frac{d^2 x}{dt^2} + \omega_0^2 x = 0, \quad \text{where } \omega_0^2 = \frac{k}{m} \quad (2)$$

We have studied the solution of this equation in the previous experiment (Exp3), where the simple pendulum underwent simple harmonic oscillation with angular frequency $\omega = \sqrt{g/L}$, where g is the gravitational constant and L is the effective length. We found that the solution for the simple pendulum was expressed in a sinusoidal form, $\Theta(t) = \Theta_0 \sin(\omega t)$. Therefore, we anticipate that the solution for Eq. (2) should be similar, namely

$$x(t) = A \cos(\omega_0 t + \varphi), \quad (3)$$

where A is the amplitude of oscillation, φ is a phase constant, and ω_0 is the angular frequency:

$$\omega_0 = \sqrt{\frac{k}{m}}. \quad (4)$$

2.2 Damped Harmonic Oscillation (DHO) The amplitude of oscillation of the mass gradually decreases over time. This is due to the effect of friction or a drag (resistive) force. We want to see if we can understand its effects. For our example, the effect of friction can be represented as a force proportional to the velocity of the mass. Therefore, Eq. (1) can be modified to read

$$m \frac{d^2 x}{dt^2} = -kx - \beta \frac{dx}{dt} \quad (5)$$

where β is a constant of proportionality, called the damping coefficient. The minus sign indicates that the damping force is always opposite to the direction of motion. Rearranging the above equation yields

$$\frac{d^2 x}{dt^2} + 2\gamma \frac{dx}{dt} + \omega_0^2 x = 0, \quad (6)$$

where $\gamma = \frac{\beta}{2m}$. The solution to this (second order differential) equation is no longer SHO. If

β is not too large, it is a modification of SHO. Also, the frequency of oscillation will be modified by the damping. The solution may be made plausible by an educated guess. We postulate that the

oscillatory motion is still sinusoidal, but that the amplitude of oscillatory motion decreases as a function of time. So our trial solution is

$$x(t) = A e^{-\gamma t} \cos(\omega t + \varphi), \quad (7)$$

where A is the amplitude, φ is the phase angle, and it turns out that $\omega = \sqrt{\omega_0^2 - \gamma^2}$ is the angular frequency of this system. See Appendix A for a more complete derivation. The coordinate $x(t)$, as a function of time t , is shown in Fig. 2.

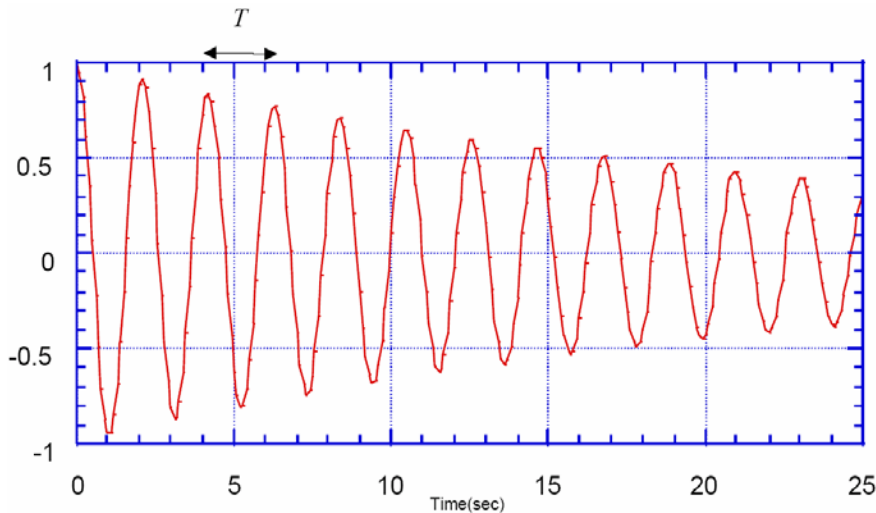


Fig. 2. Damped harmonic oscillator as a function of time. The envelope decay function is $\exp(-\gamma t)$. The period T is related to ω_0 by $T = 2\pi / \omega_0$, where $\omega_0 = 2\pi f$.

3. Questions for the Preliminary Discussion

- 3.1 In Fig. 2, what is the amplitude of the oscillation at sec $t = 10 \text{ sec}$? at sec $t = 0 \text{ sec}$?
- 3.2 Determine the period T , the angular frequency ω and the phase angle φ in $\cos(\omega t + \varphi)$.
- 3.3 If the phase angle $\varphi = \pi/2$, redraw the oscillations in Fig. 2.

Experimental Procedure

4.1 The Dynamic Force Transducer In this experiment we will use a Dynamic Force Transducer (DFT), an electronic device that outputs a voltage proportional to a force applied to it. In **Parts 8 and 9**, we will hook the outputs of the DFT to the input terminal block of a LabView card and use LabView to monitor the applied force. If we use the DFT as the spring support, then the force exerted by the spring is measured. If we ignore the mass of the spring (a constant), this force is the same as that applied to the hanging mass and is by Hooke's Law proportional to the displacement of the mass from equilibrium. Thus, by monitoring the force in the spring we are, in effect, monitoring the position of the mass. The transducer produces a voltage proportional to force, with a scale controlled by the gain knob, so the scale is unknown unless we calibrate it. Further, we use a knob to zero ("offset") the force readout for the mass of the support hook, but we use different masses during the course of the experiment, and the zeroing may be inexact. Thus, we may regard the **DFT readout as giving not x , but $Bx + C$** , and we should **add an arbitrary constant to any Kgraph equations used to fit DFT output**.

4.2 Inertial Mass. Throughout the lab, whenever measuring mass, **be sure to include the support hook** as well. (We will ignore the mass of the spring.). The key reason is that the mass in Eq. (1) is the whole mass that is accelerating—the inertial mass of $F=ma$. That of course is the same mass as enters in Eq 4 for the angular frequency.

5. Calibration of the DFT. There are two knobs on the DFT, offset (zero adjust) and gain (sensitivity). Set the gain knob nearly fully clockwise and *do not touch it again during the experiment*. Measure the hook mass, then attach the spring and hook to the DFT. Connect the output of the DFT to the digital oscilloscope. With the spring and hook mass at equilibrium, adjust the offset button of the DFT until the output reads zero volt(s) on the digital scope display. Note that zeroing the output at equilibrium is for convenience only, but it is your measurement for the hook alone, with no mass added.

1. Measure the output (volts) for three different **masses added**. The recommended order of masses is 50g, 100g, and 150g. Appendix C shows you how to measure the voltage with the oscilloscope. Assuming Hooke's law to be true (and the DFT to behave as advertised), write down the expected equation relating output to the mass added.
2. For these four measurements make a graph of V (volts) vs. m (grams). What function describes the dependence of V on m ? Is it roughly consistent with what you predicted? Do you think the DFT will do a good job of measuring position?

6. Static Measurement of the Force Constant. In this part you check Hooke's Law by measuring the displacement vs. mass using a metric tape measure or a meter stick.

1. Measure displacements for three (3) different masses. Assign uncertainties to the measured displacements and masses. Don't forget to include the mass of the hook.
2. Plot displacements (in units of cm) vs. mass (in units of grams) using *K-graph*.
3. Apply a least-squares fit procedure to obtain the spring constant k and its uncertainty. What are the units of k ?
4. Predict the oscillation frequency ω_0 for each of these masses. Show the calculation in your notebook.

7. Measurement of the Period of SHO. In this part of the experiment you will use the hand timer to measure the angular frequency of oscillation ω_0 . You will perform this for three different masses.

1. Attach your first mass to the spring and set it into oscillation. Pick an appropriate point of reference and run your timer for ten complete cycles of the oscillations. Record the time (in seconds) in your notebook.
2. Repeat step 1 for your second and third mass.
3. To obtain ω_0 , note that $\omega_0 = 2\pi f$, where f is measured in cycles per second.
4. Compare these ω_0 with those you predicted in Part 6.
5. Using *K-graph*, plot your data of ω_0 vs. m on both a logarithmic plot and a linear plot (See Appendix B on guidance on choosing which axes to make logarithmic). Do you observe the functional dependence predicted by Eq. (4)?

8. Oscillating Spring: Displacement as a Function of Time.

We will next consider the case of a mass attached to a spring as shown in Fig. 1.

1. Attach the spring-hook-mass (100g added) combination to the DFT. With the help of the digital scope, adjust the DFT offset. Set the spring-mass system into oscillation, observe the pattern on the scope display and sketch it in your notebook.
2. Measure the frequency, f , and the amplitude, A , on the digital scope using the cursors. Calculate the angular frequency ω_0 from f and record it in your notebook.
3. Disconnect the output of the DFT from the scope and connect it to the data acquisition card on the PC. Open a program called "Force Transducer.vi" in the LabView folder in C:\LabView\Vi.lib\vi_for_phy191\ForceTransducer.vi. The following specifications are needed in order to collect the data successfully. Set **device** = 1; **channels** = 0; **voltage data**: scan# = 0, channel = 1; **scan numbers** = 1000, **scan rate** = 100; **input limits** = 2. The **scan numbers** is just the number of data points; the **scan rate** is the number of points/second, so the number of seconds for which you will record data is just **scan numbers/ scan rate**. If you need a longer data sample, raise the scan numbers; if you need finer sampling, raise the scan rate.

Set the spring-mass system into oscillation and press **Ctrl R** to begin collecting the data. (The scanning process will take 10 seconds.) When it is done, the program will ask you to save the data ; save it under the folder Phy191 in your designated section folder. (Or you could save it on your floppy disk.) **Caution:** When exiting the program, it will ask you if you want to save the current setting. Choose **NO**.

4. Use *K-graph* to retrieve the data file. When inputting your data to *K-graph*, the following specifications are required: **Delimiter** = Tab; **Number** = 1, **Line Skipped** = 0, **Options/Read Title** (No check). Your file contains two columns: the first column is the time, the second column is the output voltage. **The data is big—don' print data tables.**
5. Plot voltage vs. time and use the general curve-fit editor to perform a non-linear least squares fit procedure using the function given in Eq. (3). The corresponding *K-graph* function is $m_1 + m_2 \cos(m_3 \cdot m_0 + m_4)$; m_3 and m_4 are in radians (m_1 is the unknown voltage for $x=0$).

Notes:

- This non-linear least-squares fit function will not converge unless your initial parameters ($m_1 \dots m_4$) are realistic—especially for the angular frequency, for which you must find a way to get a good estimate. Show in your notebook how you estimated the parameters. How do the values compare to those output by *K-graph*?
- The fit is subject to the ambiguities inherent in trig functions. The parameters of the fit should be reduced to standard form using trig identities. You should wind up with a positive amplitude, a positive ω , and a phase angle in the range $[-\pi, \pi]$. A simple procedure is to first correct the sign of the amplitude (if necessary) by $-\cos(z) = \cos(z+\pi)$, then correct the sign of ω if necessary by $\cos(-z) = \cos(z)$, and finally correct ϕ by $\cos(z) = \cos(z \pm 2N\pi)$, where N is any integer. Show the steps in your notebook, and report the final values in your report or hand-written on your plot (which should show the original fit values and uncertainties)

9. The Effect of Friction

In this part of the experiment, we will investigate the behavior of the spring oscillations under the effect of friction.

1. Weigh the friction umbrella. Attach the friction umbrella to your spring-hook-mass (~150g) combination and once again observe the spring tension on the digital scope. Adjust the **offset** on the DFT if necessary. Set the system into oscillation and observe the oscillation pattern on your digital scope. You should observe a pattern similar to that in Fig. 2, which shows the sinusoidal oscillation modulation by the exponential damping factor.
2. Measure the frequency f on the digital scope using the built in buttons. Calculate the angular frequency ω and the period T from f and record these values in your notebook.
3. With the friction umbrella attached to your spring-mass system as in step 1, repeat steps 3 and 4 of Part 8. Fit the data to the function given by Eq. (7): $m_1 + m_2 \exp(-m_3 \cdot m_0) \cos(m_4 \cdot m_0 + m_5)$. Show how your $m_1 \dots m_5$ parameters are estimated. How do these values compare to those given by *K-graph*? As above, standardize the amplitude, angular frequency, and phase.
4. Repeat the experiment with a 100g mass added.

10. Questions to be Discussed

1. Define and explain briefly the meaning of the terms (a) restoring force, (b) free oscillation, (c) simple harmonic motion, (d) phase angle and (e) natural frequency.
2. In Part 7 we calculated the spring constant k from the slope of a plot of displacement vs. mass. Using the same data, give an alternative method for making a plot that yields the spring constant k directly from the slope. Hint: What kind of plot would have a slope that is the spring constant k ?
3. In this question, you are asked to estimate the decay factor directly from a plot of voltage vs. time. Draw a smooth curve connecting the decay peaks on your plot in Part 9. The envelope decays as $\exp(-\gamma t)$. When $t = 1/\gamma$ the amplitude $x(t)$ has decreased by a factor of $1/e$, or about 0.368 times the initial value A . Beginning from any point on the time-axis, determine the length of time required for the amplitude $x(t)$ to decrease by $1/e$ i.e. $1/\gamma$. Compare γ with that given by *K-graph fit*.
4. Compare γ obtained Part 9 for different masses. Does γ depend on the mass of the system? Calculate γ supposing you had added 200g, and compare it to the values for mass of 100g and 150g.
5. Extra Credit: Compare the measured angular frequency of free oscillation ω_0 (found in Part 7 or 8) with the damped oscillator frequency ω (found in part 9). Again use data with added $m = 100g$. Discuss your comparison of ω for the damped oscillator with ω_0 for free oscillator. Is the effect of the mass of the friction umbrella on ω_0 significant in this comparison? Are your results roughly consistent with Eq A5 for ω and the value of γ from your fit?
6. Extra Credit: What is the % difference between your angular frequency of free oscillation ω_0 and your prediction in part 6? A more detailed theory suggests adding $M_{\text{spring}}/3$ to the suspended mass when calculating the frequency. The spring weights about 9 grams. Does this fit your observed ω_0 better than just using the suspended mass? Are either prediction compatible with your measurement with uncertainties?
6. Summarize your various measurements in two tables, one for angular frequencies ω and ω_0 , and another (much shorter) for estimates of γ . The ω table should include the method, mass, value, prediction, and difference in %; base the predictions on the part 6 results.

Appendix A: Solution of Damped Harmonic Oscillator

The solution Eq. (7) for the damped harmonic oscillator in Eq. (6) can be found as follows.

$$\frac{d^2x}{dt^2} + 2\gamma \frac{dx}{dt} + \omega_0^2 x = 0. \quad (\text{A1})$$

The general solution is

$$x(t) = Ae^{Qt}, \quad (\text{A2})$$

where Q is a factor to be determined. Then

$$\frac{d}{dt}x(t) = QAe^{Qt} = Q \cdot x(t) \quad \text{and} \quad \frac{d^2}{dt^2}x(t) = Q^2 Ae^{Qt} = Q^2 \cdot x(t)$$

Substitution of these two in (A1) yields

$$(Q^2 + 2\gamma Q + \omega_0^2)x(t) = 0, \quad (\text{A3})$$

which has a solution of the form

$$Q = \frac{-2\gamma \pm \sqrt{4\gamma^2 - 4\omega_0^2}}{2} = -\gamma \pm \sqrt{\gamma^2 - \omega_0^2}. \quad (\text{A4})$$

For our case, which involves a weakly damped oscillation, $\omega_0 \gg \gamma$. Therefore, our solution (A4) becomes

$$Q = -\gamma \pm i\sqrt{\omega_0^2 - \gamma^2} = -\gamma \pm i\omega \quad \text{with} \quad \omega = \sqrt{\omega_0^2 - \gamma^2} \quad \text{and} \quad \omega_0^2 = k/m \quad (\text{A5})$$

Our general solution Eq. (A2) now becomes

$$x(t) = A_1 e^{-\gamma t + i\omega t} + A_2 e^{-\gamma t - i\omega t}. \quad (\text{A6})$$

Set $A_1 = \frac{1}{2} A e^{i\varphi}$ and $A_2 = \frac{1}{2} A e^{-i\varphi}$, where φ is a phase factor. Then

$$x(t) = \frac{A}{2} e^{-\gamma t} e^{i(\omega t + \varphi)} + \frac{A}{2} e^{-\gamma t} e^{-i(\omega t + \varphi)} = A e^{-\gamma t} \frac{1}{2} (e^{i(\omega t + \varphi)} + e^{-i(\omega t + \varphi)}). \quad (\text{A7})$$

But $\cos(\theta) = \frac{e^{i\theta} + e^{-i\theta}}{2}$. Thus,

$$x(t) = A e^{-\gamma t} \frac{1}{2} (e^{i(\omega t + \varphi)} + e^{-i(\omega t + \varphi)}) = A e^{-\gamma t} \cos(\omega t + \varphi), \quad (\text{A8})$$

which is the solution Eq. (7).

Appendix B: Plot Axis Selection and Guessing Functional Relationships

In Kgraph, from Plot | Axis Options you can choose Linear (default) or Log (logarithmic) for either the x or y axis. Let's see how this might help you: For this discussion, we will assume logarithmic means natural logarithms, log to the base e.

1) Semi-log plots

Suppose the true relationship is something like $y = y(x) = e^{-ax}$ (the constant a might be positive or negative). If you plot y on a log scale and x on a linear scale, this is a semi-log plot (logarithmic on only one axis). Then you would be looking at $\text{Ln}(y)$ vs. x , and substituting the true relationship $y(x)$, your plot would display $-ax$ vs. x .

That is, **exponential relationships produce straight lines on a semi-log plot** and the slope is equal to the coefficient of x in the exponent

To find the slope, you could perform either the nonlinear fit $y = \text{Exp}(-a x)$, or a linear fit to the semi-log plot, fitting $\text{Ln}(y) = -a x$. These fits will give somewhat different answers as equal weighting in y assumes different uncertainties than equal weighting in $\text{Ln}(y)$.

2) Log-Log plots

Suppose now the relationship between y and x is of the form $y^p = a x^q$. A linear relationship ($y = a x$, with no additive constant) is a special case with $p = q = 1$. Assuming y and x are positive, you can take the p -th root and recast of this relationship as $y = y(x) = b x^u$, where $b = a^{1/p}$ and $u = q/p$.

If you plot both x and y on a logarithmic scale, it's a log-log plot. Then you would be looking at $\text{Ln}(y)$ vs. $\text{Ln}(x)$, which, on substituting the true dependence $y(x)$ would yield $\text{Ln}(b) + u \text{Ln}(x)$ vs. $\text{Ln}(x)$

That is, **power relationships of the form $y^p = a x^q$ produce straight lines on a log-log plot.**

The slope is the power u , to which x is raised to give y .

You can find this slope by rise over run by measuring distances directly on the graph:

$$u = [\text{Ln}(y_1) - \text{Ln}(y_0)] / [\text{Ln}(x_1) - \text{Ln}(x_0)]$$

This works directly from the plot, as long as the y and x scales are the same. It's the same as

$$u = \text{Ln}(y_1/y_0) / \text{Ln}(x_1/x_0)$$

which makes it obvious that the units cancel, and all that matters is that a given factor (e , or 10, say) takes up the same amount of space on the y and x axes. You could also use either formula with coordinates of points on the line.

Or you could perform either the nonlinear fit $y = b x^u$, or a linear fit to the log-log plot, fitting $\text{Ln}(y) = b + u \text{Ln}(x)$. These fits will give somewhat different answers as equal weighting in x or y assumes different uncertainties than equal weighting in $\text{Ln}(x)$ or $\text{Ln}(y)$.

To really test these relationships well, you should ideally:

- Compare the fit quality with an alternative, such as a straight line

- Have more than 3 data points

- Have data covering at least a factor of e^2 (≈ 7.4) to 10 in range of x and y .

Appendix C: Oscilloscope hints

An oscilloscope (scope) is a device to perform real-time visualization of voltages in an electric circuit. The screen of the scope shows the dependence of the measured voltage on time. If the voltage changes, the picture on the screen changes as well. The time is shown on the x-axis, the voltage on the y-axis. Just as in *KGraph*, you can choose the best presentation of the graph by adjusting various scope controls. And just with *Kgraph* presentation settings, **no matter what scope settings you change, you have not changed anything about what happens in the circuit.** If you know how many seconds (volts) correspond to one tick on the grid of the screen, you can read the period and the magnitude of the signal directly off the screen, most accurately by using the cursor buttons.

Setup of the oscilloscope

Your scope can simultaneously measure two voltages; however in this lab we will study only one input signal (supplied through channel A). Both channels A and B have similar sets of controls. During the lab, pay attention that you only work with the controls for channel A, not channel B. The current settings for channel A (B) are shown in the upper part of the digital display to the right of the screen. By pressing the "A/B" button make sure that the shown settings correspond to channel A. The only settings you should have to change on the scope are those on the top row, the TB setting at the middle towards the left, and the row second from the bottom.

Connect a BNC cable from the transducer to the channel A input. Turn on the scope. Press "AC/DC" button until the digital display shows "DC". Also, turn the "DIGITAL MEMORY" on. Choose the time base "TB" (the value of the 1 cm grid of the time scale) to be equal to "0.5 s", and the voltage base "V – mV" to be equal to "0.1 V". During the lab, you may want to adjust these settings to get the best presentation of the graph. It's also nice to push XMAG 6 times or so until the display shows [-----] instead of [----] so the scope doesn't wait to refresh the data on the screen.

Press "GND". After this, the scope input will be connected to ground, so that the input voltage is 0 V. Turn the Y-pos and X-pos knobs until you position the green line of the signal exactly in the middle of the screen. After that, do not touch the X- pos and Y-pos buttons for the rest of your measurements. This way, you will ensure the correct calibration of the picture offsets inside the scope (then you will only need to worry about the voltage offset inside the force transducer which you will adjust with the help of "ZERO" on the box of the transducer).

Press "GND" again. Now you will see the signal coming from the transducer. For nonzero mass, it will be vertically shifted with respect to the center of the screen. For static stretching of the spring, you will measure this displacement of the signal to find the dependence between the applied mass, elongation of the string and the voltage. During the study of oscillations, you will have to adjust the "ZERO" on the force transducer until the signal line is again in the middle. Thus you will make sure that the shown voltage corresponds to the displacement of the mass from the equilibrium point, and not from some other point of the system.

Several other useful buttons: you may use the knobs to the left of the screen to adjust the focus and the brightness of the line. The button "LOCK" will allow you to lock (freeze) the image, which makes detailed measurements on the captured image convenient, as described below.

Voltage and Time Measurements with the Cursors

Inspect the main screen. If there is no writing at the bottom of the screen, press one of the blue "soft-keys below the screen to make the writing appear. If there is some writing and one of the soft-keys has RETURN written above it, keep pressing it until it no longer says RETURN, to get back to the top menu. You should see:

CURSORS

SETTINGS

TEXT OFF

Reading Voltages. We will use these soft-keys to move screen cursors so they correspond to the size of our signal and read off voltages from the screen. Press the CURSORS soft-key. Press MODE to set up the cursors we want. Toggle the V-CURS and T-CURS soft-keys until the horizontal cursor lines are on and the vertical cursor lines are off. Press RETURN.

Press V-CTL to control the Voltage, cursors. The V cursors measure voltage, which is displayed as the vertical (y) axis on the screen. They appear as horizontal lines on the screen, representing a constant voltage level. Move the top one down to the level of your signal, and move the bottom one up reference value. The value is displayed at the top of the screen. For uncertainty estimation: use 2 clicks, or two linewidths, or remeasure a few times.

Frequency Measurement: Next, we wish to determine the time between 2 points. Time is measured as a horizontal distance on the screen. The sine-wave signal is cyclic: it repeats itself. The time a sine-wave signal takes to make one cycle is called the period of the signal. It has units of seconds. The inverse of the period is the frequency (in Hz).

To measure time we need the T-cursors (T stands for time). Time is displayed as the x -axis on the screen, so a constant time is marked by the T cursor as a vertical line. Return to the top-level menu using the RETURN soft-key. Press MODE and toggle the T-CURSOR to on. Hit RETURN, then hit T-CTL to move the time-cursors to the left and right. The length of the cycle can be measured by positioning the T-cursors on equivalent points on the trace, for example, on two adjacent peaks or two adjacent valleys. The period and the frequency can both be read off from the top of the screen.

